

wellogic

RoutinePair Skincare Device

5 WAYS TO ELEVATE YOUR SKINCARE ROUTINE




Clean

Cleansing Mode: Unlocks a deeper pore cleanse right after facial cleansing

Apply toner or essence to the cotton pad, then secure it in place using the retaining ring.

Boost

Mask Mode: Elevates your skincare experience while using a sheet mask

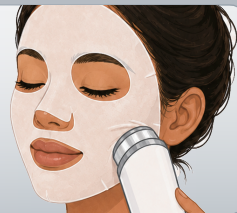


Lifting

Lifting Mode: Delivers firming and contouring care after applying serum

Eye Care

Eye Care Mode: Dedicated gentle care for the eye area after applying eye cream

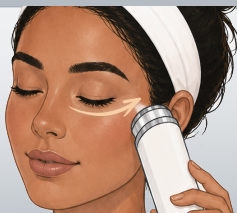

Cool

Refining Mode: Enjoy a soothing cooling finish to wrap up your skincare routine

MATCH MODE TO YOUR SKIN STEP



CLEAN

Complete your cleanse.

Daily, as needed



BOOST

Help skincare products absorb more effectively.

Morning and evening



LIFTING

Support a firmer, more defined look.

2-3 times per week



EYE CARE

Refresh the delicate eye area.

Morning and evening



COOL

Cool and soothe your skin.

As needed

Gentle modes daily; Lifting 2-3 times weekly


Before you start



Start with clean skin.



Do not use on broken skin, open wounds, or irritated areas.

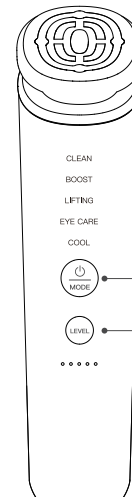

1. Press and hold the "MODE" button for 1.5 seconds to power on.
2. Short press the "MODE" button to select a mode. Each short press will switch Wellogic RoutinePair Device through the following modes in a continuous cycle: CLEAN-BOOST-LIFTING-EYE CARE-COOL.
3. Short press the "LEVEL" button to cycle through intensity levels from Level 1 to Level 5.

4 minutes per session · Auto shut-off
See the User Manual for complete safety information.

Scan QR for video guide

112x202mm